


**MIP**

LAUNDRY CART Assembly Instructions


Thank you for your purchase of the Med-i-Cart.
Before assembly, we suggest you read the following instructions.

CART ASSEMBLY

Parts included for Cart Assembly:

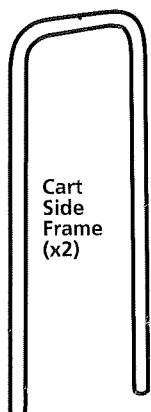

Cart Side Frame (x2)

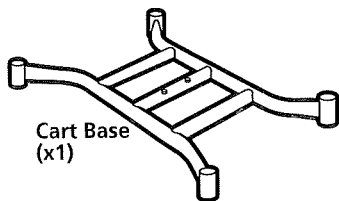

Cart Base (x1)



Accessory Bar (x2)



Bumper (x4)



Wheel without Brake (x2)



Wheel with Brake (x2)



Bag Frame

Style: LCC-2000-1-SP (x1)

Style: LCC-2000-2-SP (x2)

Style: LCC-2000-3-SP (x3)

Style: LCC-2000-4-SP (x4)



Plastic Lid

Style: LCC-2000-1-SP (x1) White


Style: LCC-2000-2-SP (x2) White, Red

Style: LCC-2000-3-SP (x3) White, Red, Blue


Style: LCC-2000-4-SP (x4) White, Red, Blue, Green

1

1. Remove all packaging
2. Use the empty box as protective mat during assembly
3. Slide the side frames into the cart base, ensuring the holes to install the accessory bars are facing inward
4. Insert the accessory bars by gently pulling the side frames, one end at a time

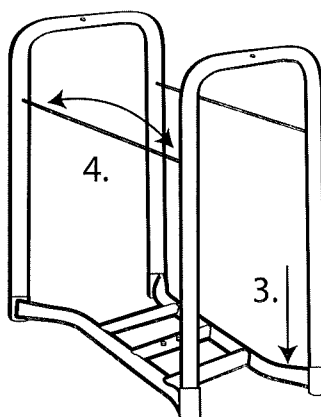**2**

1. Turn the frame upside down, ensuring that the side frames are fully inserted
2. Position the bumpers on the base tube so that the holes align
3. Insert the wheel pin through the bumper, into the side frame, push firmly or tap gently with a hammer to secure

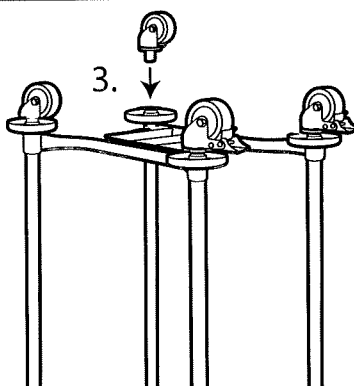**3**

1. Turn the frame the right side up
2. To insert the bag frame(s), use the holes at the top of the side frame. Align the pins under the fitting lip and snap in place
3. For carts with more than one bag frame, ensure all the lid stoppers face the same direction
4. Snap on the plastic lid(s). The handle must face the direction opposite from the lid stopper, so when the lid is raised, it cannot go past a 90 degree angle

